





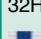
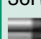







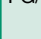




Hartmetall-Bohrer für die Zahnmedizin
 Carbide burs for the dentistry
 Fraises carbure pour le odontologie
 Fresas de carburo para odontologia

Anwendung / Application / Application / Empleo

Einsatzgebiet Used on Domaine d'applicat. Usos en sector	Arbeitsgang Working process Phase de travail Fase de trabajo	Verzahnungsart Type of cuts Type de denture Tipos de dentado	Instrumentenbeispiel Examples of instr. Exempl. d'instruments Ej. de los instrumentos	ISO-Ø	empf. Drehzahl min ⁻¹ Recom. speed r.p.m. Vit. de rotation rec. Veloc. recomendada
Schmelz, Dentin Enamel, Dentine Email, Dentine Esmalte, Dentina	Exkavieren Excavating Excavation Excavación	Schnittfreudige Verzahnung mit/ohne Spezialhieb Power cut with/without special cut Denture particulièrement sécante avec/sans taille spéciale, Dentado de alto rendimiento, con/sin talla special	1SXM / W / RA / CA  1SX / W / RA / CA  1S / W / RA / CA 	008-012 014-016 018-027	3.000 - 2.500 2.000 - 1.500 1.000 - 500
Schmelz, Dentin Enamel, Dentine	Präparieren Preparation	Mittlere Verzahnung Medium cut Denture moyenne Dentado medio	21R / FG 	005-009 010-018	225.000 - 165.000 215.000 - 62.000
Email, Dentine Esmalte, Dentina	Préparation Preparación	Mittlere Verzahnung mit Standard-Querhieb Medium cut with standard-cross-cut Denture moyenne avec taille transversale standarde Dentado medio con corte transversal	31R / FG 	008-012 014-016	225.000 - 120.000 180.000 - 80.000
Schmelz, Dentin, Füllungsmaterialien Enamel, Dentine Filling materials Email, Dentine matériaux d'obturation Esmalte, Dentina, mat. de obturaciones	Schnelleres Präparieren Faster preparation Préparation plus rapide Preparaciones rápidas	Mittel-feine Kreuzverzahnung Medium-fine double-cut Denture croisée moyenne-fine Dentado medio-fino cruzado	SPEEDIES, 7X / FG 	008-012 014-018	225.000 - 120.000 180.000 - 80.000
Amalgam Amalgam Amalgame Amalgama	Amalgam-Entfernung Amalgam remova Elimination d'amalgame Eliminación de amalgama	Mittlere Spezial-Kreuzverzahnung mit überlaufender Schneide Medium special double-cut with overflow cut Denture spéciale croisée moyenne avec une lame frontale Dentado medio especial cruzado que continua en la cabeza	32HG / FG 	012	200.000 - 180.000
Metalllegierungen Metal alloys Alliages métalliques Aleaciones metalicas	Auftrennen von Kronen/Brücken Cutting of crowns and bridges Séparation des couronnes et bridges Separación de coronas y puentes	Mittlere Verzahnung mit feinem Querhieb/ Mittlere Kreuzverzahnung Medium cut with fine cross-cut/ Medium double cut Denture moyenne avec taille transversale fine/ Denture croisée moyenne Dentado media con talla transversal fina/ Dentado medio cruzado	36R / FG  21RX / FG 	010-012	200.000 - 180.000
Wurzeloberflächen Root surfaces Surface des racines Superficie de las raices	Säubern/Glätten Scaling/root-planing Nettoyage/lissage Limpieza/alisado	Polygonales Arbeitsteil (ohne scharfe Schneiden) Polygonal working part (without sharp edges) Partie trav. polygonale (pans sans tranchants) Parte de trabajo provista polygonal (sin cortantes vivos)	407/407L / W/ RA/ CA  406/406L / W/ RA/ CA 	012 010+012	12.000 - 8.000 12.000 - 8.000
Schmelz, Dentin, Komposit, Metall, Keramik Enamel, Dentine Composite, Metal Ceramic Email, Dentine, Composite, Céramique, Métal Esmalte, Dentina Composite, Cerámica, Metálicas	Finieren Finishing Finissage Acabado	Finierer-Verzahnung mittel, fein, ultra-fein Finishing cut medium, fine, super-fine Denture de fraise à finir moyenne, fine, ultra-fine Dentado de fresas para acabar medio, fino, superfino	44E / FG  44 E-F / FG  44E-UF / FG 	008-010 012-016 018-027	70.000 - 10.000 60.000 - 10.000 35.000 - 10.000
Knochen Bones Os Huesos	Chirurgische Eingriffe Surgical operations Interventions chirurgic. Intervenciones quirúrgic.	Chirurgie-Spezialverzahnungen Surgery- special cuts Dentures spéciales pour chirurgie Dentado especial para la cirugía	141 / HST/HP/PM / W/ RA/ CA  FG/ FGx-lg 	010-050 010-023	5.300 - 800 100.000 - 15.000